


DynaCAD ENTERPRISE

XP, SP and Workstation Solutions

Multiple configurations to meet different customer needs from entry level to high end academic research sites.

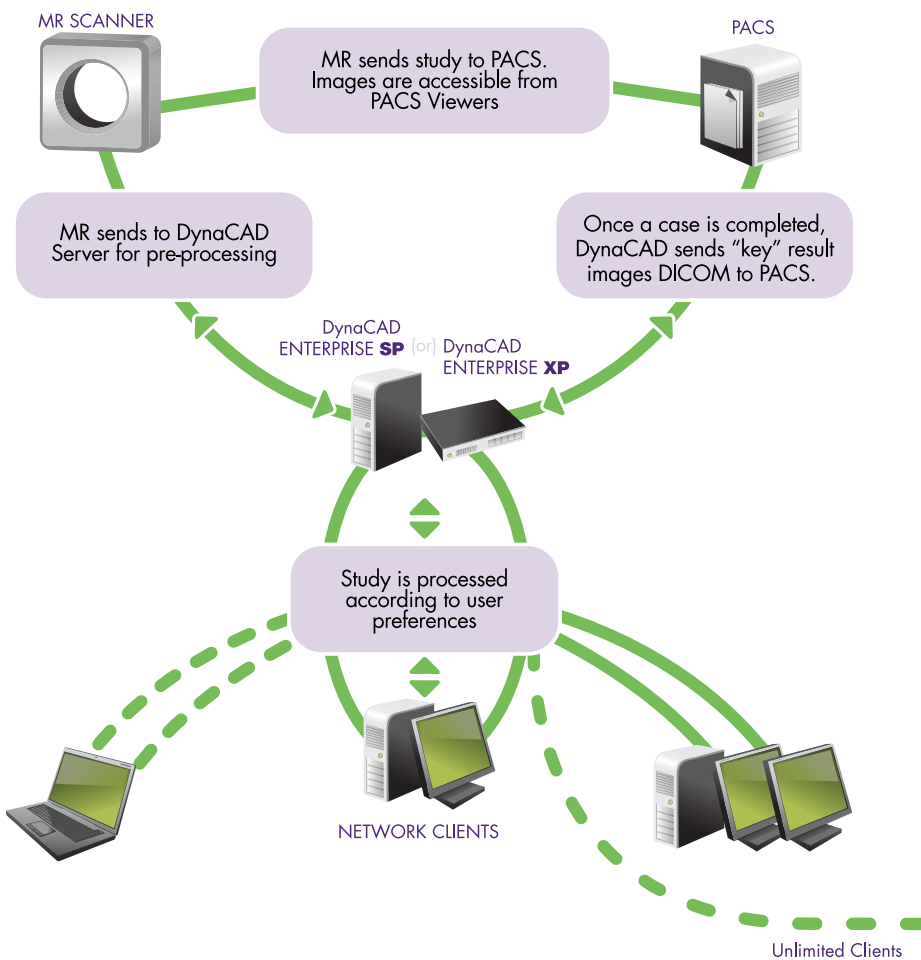

DynaCAD Enterprise XP is a high performance rack mount server configuration featuring the **DynaCAD Enterprise** Software License with EIGHT concurrent user functionality. This system has the key hardware attributes needed for the busiest departments including six RAID 5 hard drives for large data storage. It also features full redundancy of hardware components such as power supplies, hard drives, fans and more.

DynaCAD Enterprise SP is a robust client/ server implementation for the smaller volume institution that needs from 3-5 concurrent user capability within the department. For departments with workstation only requirements the standard workstation **DynaCAD** configuration is optimally configured for speed, performance and storage capacity.

Up to 8 Clients Reviewing Simultaneously

Each client independently configured for optimal reviewing performance

Connectivity through PACS viewers and other 3rd party workstation

The Enterprise Solution offers multiple, uniquely configurable connectivity options for each of the clients.

Shared Workspace with Prefetching or On Demand

With Shared Workspace and Prefetching, the patient files are sent from the MR to the **DynaCAD** Server, advanced processed, then automatically sent to the desired clients. This is ideal for the clients used to read most of the **DynaCAD** studies. Reviewing performance is similar to a stand-alone workstation.

Data Streaming with Full, Medium, & Low Settings


For clients where bandwidth or network speed are not optimal, Invivo offers Data Streaming with adjustable settings of performance. This solution is ideal for a remote office or a physician's home.

CUSTOMER SUPPLIED CLIENT PC CONFIGURATION

Recommended

- Windows XP, with Service Pack 2
- 2 GB or more RAM
- 30 GB or more of free disk space
- 128 MB video card with OPEN GL 2.0 support (nVIDIA preferred)
- Dual core 2.4 GHz processor or higher
- 100 Mbps network card
- 1680 x 1050 Monitor Resolution

DynaCAD WORKSTATION SOLUTION

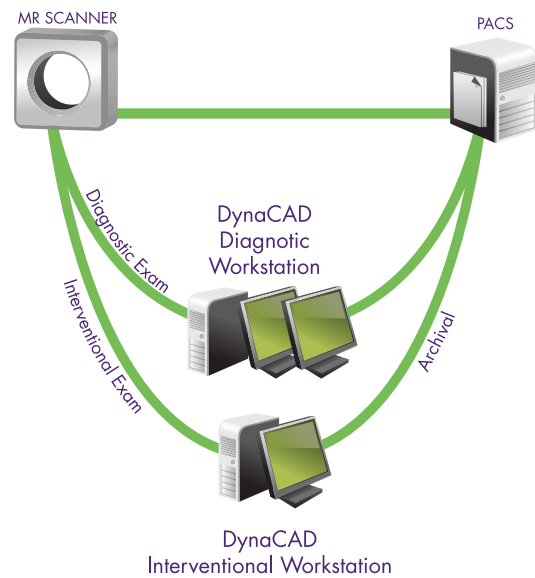

*For smaller facilities with fewer client reading station requirements, Invivo offers cost effective stand-alone workstations. As the department needs grow, your stand-alone **DynaCAD** can easily be upgraded to the **Enterprise Solution**.*



inspired by a better way

For more information about **DynaCAD** or any of the Complete Solution products from Invivo, please call us at 1-800-331-3220, or visit our website at www.invivocorp.com

Invivo Corporation

Orlando, FL U.S.A
T: 800-331-3220
T: 407-275-3220
E: info@invivocorp.com

Europe:
T: +31 (0)499 378299
contact@invivocorp.com

www.invivocorp.com

*Products, services, and technical specifications are subject to change without notice.
©2009 Invivo Corporation All Rights Reserved. 580230 Rev 01*